



PEARL

Setup Guide

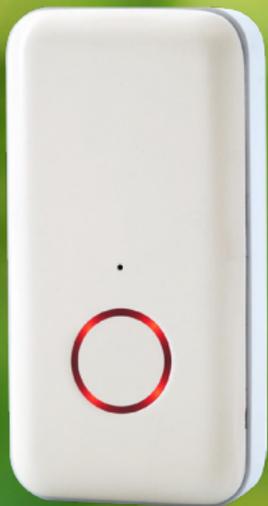

LifeStation
MEDICAL ALERT SYSTEMS

What's in the box

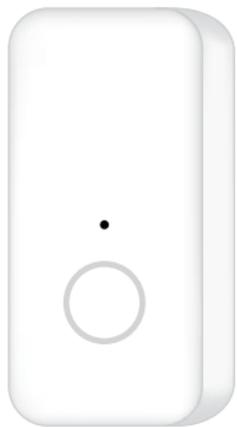

Pearl



Charger




Lanyard



Belt Clip
(optional)

Get to know your **Pearl**

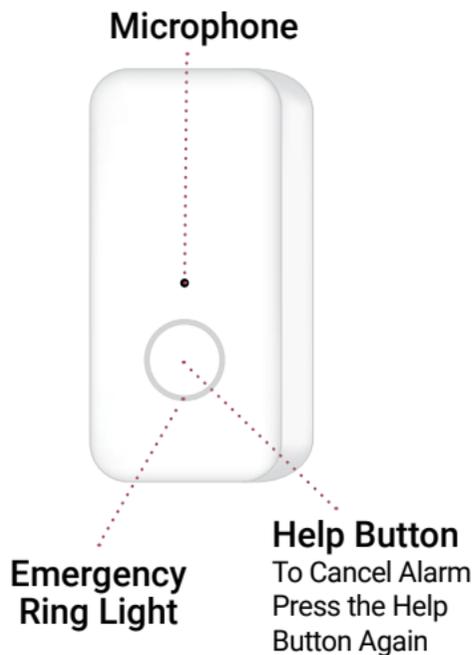


Get to know your Pearl

Checking your battery level:

Press the Power Button on the right side of your Pearl.

The green lights will indicate your devices remaining battery level.

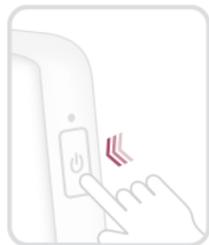

1 green light

**CHARGE
IMMEDIATELY**



2 green lights

**PARTIALLY
CHARGED**



3 green lights

**FULLY
CHARGED**

Frequently asked questions

1. How often should I charge my Pearl?

For best results, you should fully charge your Pearl daily.



Charge Daily

2. Is my Pearl waterproof?

Your Pearl is IP67 waterproof and can be worn in the shower but should not be submerged in water.



**Water-
Resistant**



**Do Not
Submerge**

3. Where does my Pearl work?

Your Pearl will work both inside and outside your home, so long as there is cellular coverage.



**In the
House**



**Out of the
House**

4. How often should I test my Pearl?

You should test your Pearl monthly. To test your Pearl, press and hold the Help Button until you hear the emergency alert message. When you are connected to the Monitoring Center, indicate that you are just testing your Pearl.



**Test
Monthly**

Setting up your Pearl

1. Plug in the charger.

The charger light will turn white.



Charger plugged in

2. Place your Pearl on the charger as shown.

When charging properly, the charger light will turn **blue** and your Pearl will announce "Charging".



Pearl charging

The battery status lights will turn **green** while it is charging. Your Pearl is fully charged when all three **green** lights are visible. Your Pearl may take up to four hours to fully charge.



Pearl fully charged

3. Test your Pearl by pressing and holding the Help Button until you hear the emergency alert message. When you are connected to the Monitoring Center, indicate that you are just testing your device.

To power ON your Pearl:

Place your Pearl in the charger or press and hold the Power Button for three seconds. When your Pearl powers on, you will feel a vibration and hear "Powering on."

To power OFF your Pearl:

Press and hold the Power Button for three seconds. Your Pearl will vibrate, the power light will flash and then turn off. **Please make sure to turn off your Pearl before returning it.**

In case of emergency

If you need assistance, press and hold the Help Button until you hear the emergency alert message.

You will be connected with a Care Specialist who is available to help.
Your Pearl will work inside or outside, wherever there is cell coverage.

How to cancel an accidental activation

To cancel an alarm, press the Help Button again.

If you are connected, and the Monitoring Center asks if you need help, state in a clear voice "No."
You will be asked to confirm that no help is needed.



Fall Detection (optional)

For devices with Fall Detection enabled, if a fall is detected and you are unable to press the Help Button, the device sends an emergency fall alert to a Care Specialist.

For best results, wear your Pearl on the lanyard and outside of all your clothing.

Fall Detection does not detect 100% of falls. If able, users should always push their button when they need help.

How to wear your Pearl

Wear your Pearl as a necklace

Your Pearl will come ready to wear with the lanyard attached. Should you take the lanyard off, it can be easily reattached by looping it around the hook at the top of your Pearl. The lanyard is adjustable for your comfort and has a breakaway clasp.



Please note if you have Fall Detection, we recommend you wear your Pearl around your neck.

How to adjust your lanyard

1. Adjusting the length

Your lanyard can be made shorter or longer. It can be adjusted between 17 inches and 34 inches to fit you comfortably.



2. To make it tighter

Hold the fastening loop with your fingers and gently pull away from the endcap. This will make the lanyard shorter and help it sit snugly around your neck.

3. To make it looser

Hold the fastening loop and slowly move it towards the endcap.

This will make the lanyard longer and allow for more room around your neck. Once you have reached your desired length, pull the endcap to secure the fit.

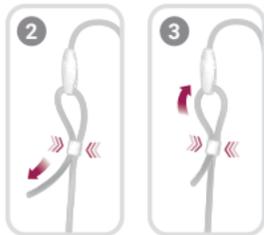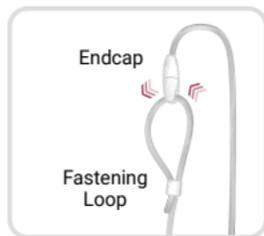

Using the optional belt clip

- If your package included the belt clip add-on, hold the belt clip in one hand and the device in your other hand. Place your Pearl into the clip and push it into place until you feel it click.



Thank you for joining the
LifeStation family!

LifeStation[®]
MEDICAL ALERT SYSTEMS

Customer Care
800-998-2400



Scan here to view this guide in a
larger format or to read in Spanish,
Mandarin, Russian, and more.

This device is dependent on cellular coverage to work.
Where cellular coverage is limited, you may encounter diminished device capabilities.